

# User manual for Bestic Picasso MINI



## Table of content:


- pp. 2-3    How to use the Bestic Picasso MINI control with the Bestic main unit
- pp. 4-7    How to use Bestic Picasso MINI control *as a computer mouse* with a Windows PC
- pp. 8-11    How to use Bestic Picasso MINI control *as a computer mouse* with a Mac computer

# How to use the Bestic Picasso MINI control with the Bestic main unit

Compatible with Bestic software from version: 1.01  
For advanced settings, set number of buttons to: 5

## Requirements

1. Bestic main unit
2. Bestic Picasso MINI
3. USB cable

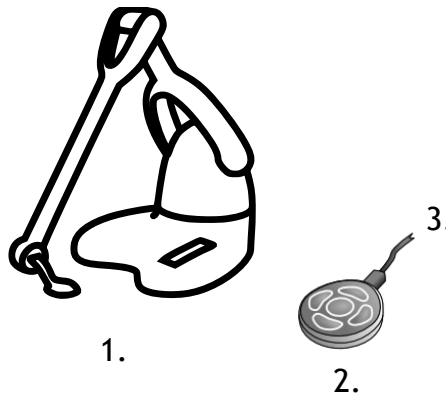

## Setup

1. Connect the USB cable to Bestic Picasso MINI as shown in the picture above. Note that the USB cable has different connector sizes on each of its two ends. The smaller one connects to Bestic Picasso MINI.
2. Turn the Bestic main unit on by pressing the black button.
3. Plug in the other end of the USB cable to Bestic.
4. While roaming the menus and eating with Bestic, use the directional arrows on Picasso mini to change the values and steer the spoon. Use the round button in the middle for "OK" and for picking up the food.

*Note: Do not use any extra button when controlling Bestic with Picasso MINI. The extra button feature is for computer use. For more information, see the user manuals for using Picasso MINI with computers.*

*You can change the USB-cable (USB A to USB micro B) for a different length cable than the standard one up to 2 m. Those cables can be purchased at computer accessory stores.*

## Safety information

Not to be used in contact with water.

Do not keep close to any fire.

Do not open.

## Cleaning instructions


To be cleaned with a moist cloth needed.

# How to use Bestic Picasso MINI control as a computer mouse with a Windows PC

## Requirements

Here's a list of what you require:

1. Bestic Picasso MINI
2. USB Cable
3. Empty USB socket on your PC
4. Extra button

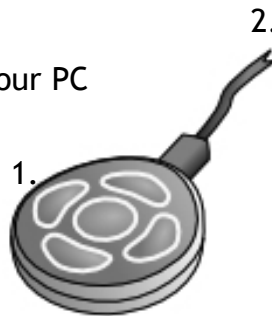


## Setup

1. Not necessary, but gives extra functions: Connect the Extra button to Picasso.
2. Connect the USB cable to Picasso as picture above shows. Notice that the USB cable has different connectors on each of its two ends. The smaller one connects to Picasso.
3. Plug in the other end of the USB cable to your empty USB socket.
4. Your computer automatically recognizes Picasso as a mouse.

## Using Picasso

### Extra button not pressed:



Press the four directional buttons to move the mouse pointer

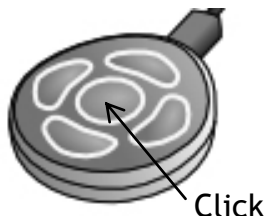

Move the mouse around to move the mouse pointer

The 4 directional buttons moves the mouse pointer (you can press up and right at the same time if you would like to move diagonal).  
The round button in the middle - click on a normal mouse.

### Extra button pressed at the same time as one of the other buttons on Bestic Picasso MINI:

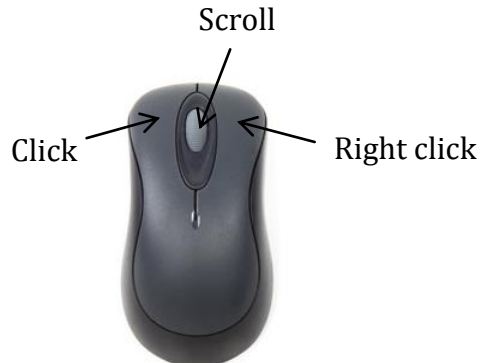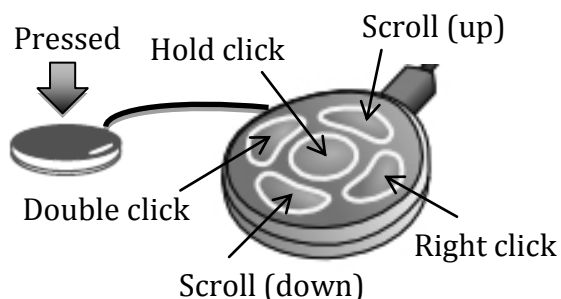

Up and down - scrolls a document or a homepage.

Left - double click on a normal mouse.

Right - Right click on a normal mouse.

The round button in the middle - click is pressed down and held down until you press the middle button again (with or without the extra button pressed).

## Tip

### Mark, copy, cut and paste

1. Put the mouse pointer in the beginning of what you want to mark and hold the extra button down at the same time as you press the round button in the middle.
2. Release the middle button and the extra button.
3. Move the mouse pointer with the four directional buttons until you have marked the wanted area and then press the round button in the middle.
4. Cut (ctrl +x) or copy (ctrl + c) the chosen items and paste it back where you want it (ctrl + v).

### Scroll

When a document or homepage is too long to fit your screen it is good to scroll down while reading.

1. Put the mouse pointer on the document and press the button in the middle.
2. Press the extra button at the same time as up or down on Picasso and the document will move up or down.

### Right click

When right clicking on a document a window pops up with alternatives of actions that you can do. Which actions that pops up depends on what you have marked and is usually very helpful. If you mark a text you can choose copy or cut for example.

1. Put the mouse pointer on the document where you want to do something and press the button in the middle to choose the document, or mark a part that you would like to work with.
2. Press the extra button at the same time as right on Picasso and a window will pop up.

### Change speed for the mouse pointer

Open the “Control Panel”.

Click on “Mouse”.

Under the tab “Pointer Options” you can choose the speed for the mouse pointer. Put the mouse pointer over the marker until its blue, press the extra button, and the round button in the middle. Move the marker to the left for higher speed and to the right for lower speed.

### Safety information

Not to be used in contact with water.

Do not keep close to any fire.

Do not open.

### Cleaning instructions

To be cleaned with a moist cloth soaked with soap, if needed.



*You can change the USB-cable (USB A to USB micro B) for a different length cable than the standard one up to 2 m. Those cables can be purchased at computer accessory stores.*

# How to use Bestic Picasso MINI control with a Mac computer

## Requirements

Here's a list of what you require:

1. Bestic Picasso MINI
2. USB Cable
3. Empty USB socket on your Mac
4. Extra button

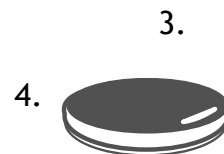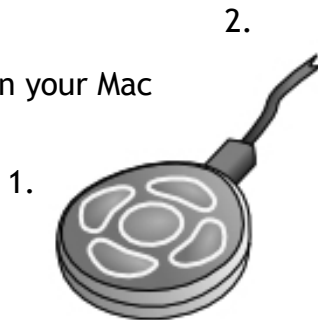

## Setup


5. Not necessary, but gives extra functions: Connect the Extra button to Picasso.
6. The first time you intend to use Picasso you should be logged in to your user on your Mac.
7. Connect the USB cable to Picasso as picture above shows. Notice that the USB cable has different connectors on each of its two ends. The smaller one connects to Picasso.
8. Plug in the other end of the USB cable to your empty USB socket.
9. Your Mac automatically recognizes Picasso as a mouse.



## Using Picasso

### Extra button not pressed:


Press the four directional buttons to move the mouse pointer

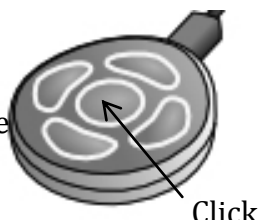

Move the mouse around to move the mouse pointer

The 4 directional buttons moves the mouse pointer (you can press up and right at the same time if you would like to move diagonal).  
The round button in the middle - click on a normal mouse.

### Extra button pressed at the same time as one of the other buttons on Bestic Picasso MINI:

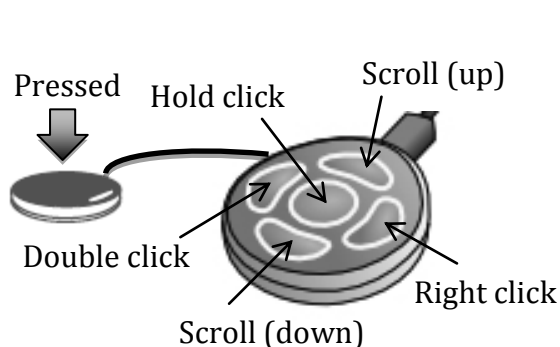

Up and down - scrolls a document or a homepage.

Left - double click on a normal mouse.

Right - Right click on a normal mouse.

The round button in the middle - click is pressed down and held down until you press the middle button again (with or without the extra button pressed).

## Tip

### Mark, copy, cut and paste

5. Put the mouse pointer in the beginning of what you want to mark and hold the extra button down at the same time as you press the round button in the middle.
6. Release the middle button and the extra button.
7. Move the mouse pointer with the four directional buttons until you have marked the wanted area and then press the round button in the middle.
8. Cut (cmd + x) or copy (cmd + c) the chosen items and paste it back where you want it (cmd + v).

### Scroll

When a document or homepage is too long to fit your screen it is good to scroll down while reading.

3. Put the mouse pointer on the document and press the button in the middle.
4. Press the extra button at the same time as up or down on Picasso and the document will move up or down.

### Right click

When right clicking on a document a window pops up with alternatives of actions that you can perform. Which actions that pops up depends on what you have marked and is usually very helpful. If you mark a text you can choose copy or cut for example.

5. Put the mouse pointer on the document where you want to do something and press the button in the middle to choose the document, or mark a part that you would like to work with.
6. Press the extra button at the same time as right on Picasso and a window will pop up.

## Safety information

Not to be used in contact with water.

Do not keep close to any fire.

Do not open.

## Cleaning instructions

To be cleaned with a moist cloth soaked with soap, if needed.

*You can change the USB-cable (USB A to USB micro B) for a different length cable than the standard one up to 2 m. Those cables can be purchased at computer accessory stores.*